


ZIRLUX® NIGHTGUARD & SPLINT – FLEXIBLE

Comfort Meets Durability

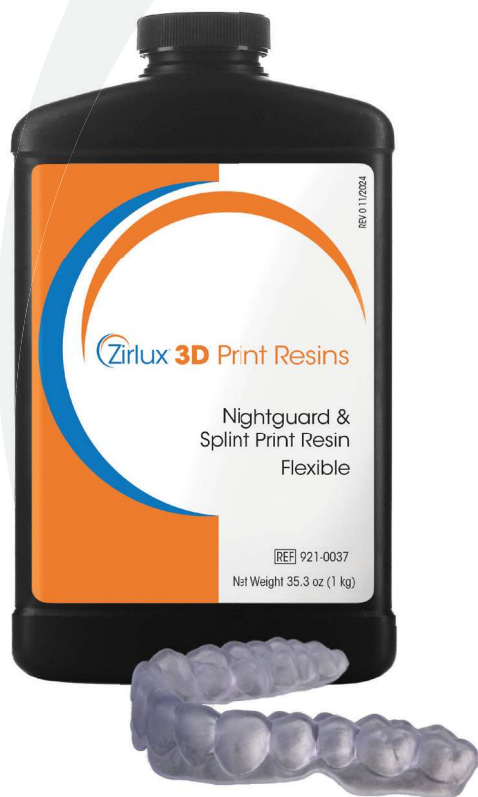

1Kg Bottle
(921-0037)

Zirlux® Nightguard & Splint Resin – Flexible is a 510K cleared printing resin that is indicated for the fabrication of nightguards, splints, mouthguards, repositioners, and bleaching trays. This DLP-optimized resin is an exceptional printing material that combines the strength needed to protect teeth against disorders like bruxism with added flexibility for better patient comfort.

KEY FEATURES:

Flexible Material: Strong but not brittle resin ideal for splints and nightguards

Durable: Stain resistant material maintains integrity over time

Biocompatible: FDA 510K cleared

Polishable: Achieves a smooth, glossy finish for patient comfort and esthetics

Extensive Printer Validation: Designed for seamless compatibility with industry-leading 3D printers.

Esthetic Finish: Smooth, clear surfaces reduce polishing time

Increased Patient Acceptance: BPA, MMA & PHTHALATE FREE

Property	Procedure	Zirlux® Nightguard & Splint Resin – Flexible
Flexural Strength	ISO 20795-2	2.6-4.4 M
Flexural Modulus	ASTM D790	1356
Tensile Strength	ASTM D638	52
Tensile Young's Modulus	ASTM D638	1790
Water Sorption	ISO 20795-2	< 18 µg/mm ³
Elongation at Break	ASTM D638	110%
Biocompatibility	ISO 10993-1 & USP Class VI	Pass



Explore the rest of
the Zirlux® 3D Print
Resins today

Zirlux® 3D Print Resins

Discover why Zirlux® is *right for you.*
www.zirlux.com

© 2025 Henry Schein, Inc. No copying without permission. Not responsible for typographical errors. 24ZS8614

Exclusively Available from

